


# AiREX Smart Glove Tester®



## Fast and Easy Operation

### What's Smart Glove Tester®?

- An advanced, fully automated glove leak testing device with pressure decay measurement.
- A faulty glove assembly presents a route of contamination. By using Airex glove leak tester, small pinholes that cannot be found through visual inspection can be detected.

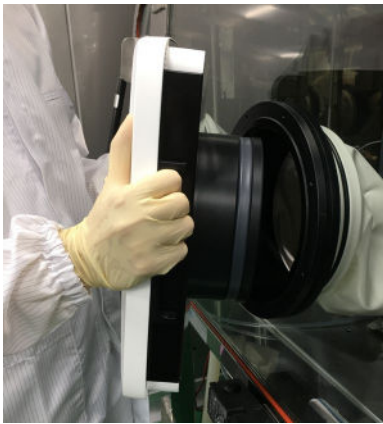

### Benefits

1. Fully automated testing
2. Easy-to-handle design
3. Customizable for any ports
4. Pressure decay testing method
5. Self-sufficient system without any external pumps, tubes, or wires
6. Wireless communication through Wi-Fi
7. 21 CFR Part 11 compliant
8. Complete validation and documentation included
9. Easy to clean
10. Fast, reliable testing time!

### Specifications

Specifications	
Utility	Power AC100V~240V 1.4A max
Battery Charger	Lithium Ion Battery and dedicated charger INPUT AC100-240V 1.4A
Dimension	(W)260mm×(D)138~183mm×(H)360mm
Operating System	Windows 10, Windows 11
Network Communication	Method: Wi-Fi Area: Max 20m inside, 30~50m outside
Simultaneous Measurement	Max 50
Port Registration	Max 99 (Duplicate registration prohibition function)
Security	4 levels 3 Operator levels, 1 Administrator level
Log and Alarm Management	Logs and alarms are encrypted and can be viewed with the dedicated PC software
Date Integrity	Encrypted files cannot be modified

CONTACT OPIRA NOW

**1300 157 969**

SOLUTIONS@OPIRA.COM.AU


OPIRA BRISBANE  
32 DIVIDEND ST  
MANSFIELD, QLD 4122

OPIRA MELBOURNE  
25 GRAHAM RD  
CLAYTON SOUTH, VIC 3169

OPIRA PERTH  
9 EARLSTON PLACE  
BOORAGOON WA 6154

OPIRA SYDNEY  
20 DUKE ST  
FORESTVILLE, NSW 2087



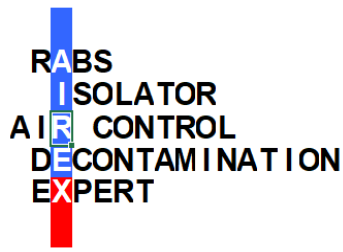

# Smart Glove Tester® WIFI Key Features

✓ **Simultaneous and parallel me:**

Max 50 gloves controlled by the SGT so allows for timesaving glove integrity te:

✓ **Fast glove integrity verification**

1 minute integrity test for each glove

✓ **Status information**

Pressure, battery status, measurement information and testing results are displayed

✓ **Wireless Communication**

WiFi communication to PC software

✓ **Wireless Communication**

WiFi communication to PC software

✓ **100% Mobility**

No additional air supply needed with automatic glove filling pump

✓ **Customizable to any port**

Customizable for different ports. Airex has experience with Bosch, GEA, Bausch and Strobel ports

✓ **Easy Operation**

All operation performed by the Smart Glove Tester® touch screen.

✓ **21 CFR Part 11**

Software Compliant with 21 CFR Part 11



CONTACT OPIRA NOW

**1300 157 969**

SOLUTIONS@OPIRA.COM.AU

OPIRA BRISBANE

32 DIVIDEND ST  
MANSFIELD, QLD 4122

OPIRA MELBOURNE

25 GRAHAM RD  
CLAYTON SOUTH, VIC 3169

OPIRA PERTH

9 EARLSTON PLACE  
BOORAGOON WA 6154

OPIRA SYDNEY

20 DUKE ST  
FORESTVILLE, NSW 2087

